

Pre-lab

The goals of this pre-lab are to:

- Gain an understanding of the SEA-AD project.
- Build confidence in accessing the SEA-AD interface and exploring this vast dataset.

What is the SEA-AD project?

First, let's get to know a bit more about the SEA-AD project, which was reported by (Gabbito et al. 2024) in this [Nature Neuroscience paper](#).

We're not going to read the whole paper, but let's take a look at part of Figure 1 (Panel A), where they outline the project:

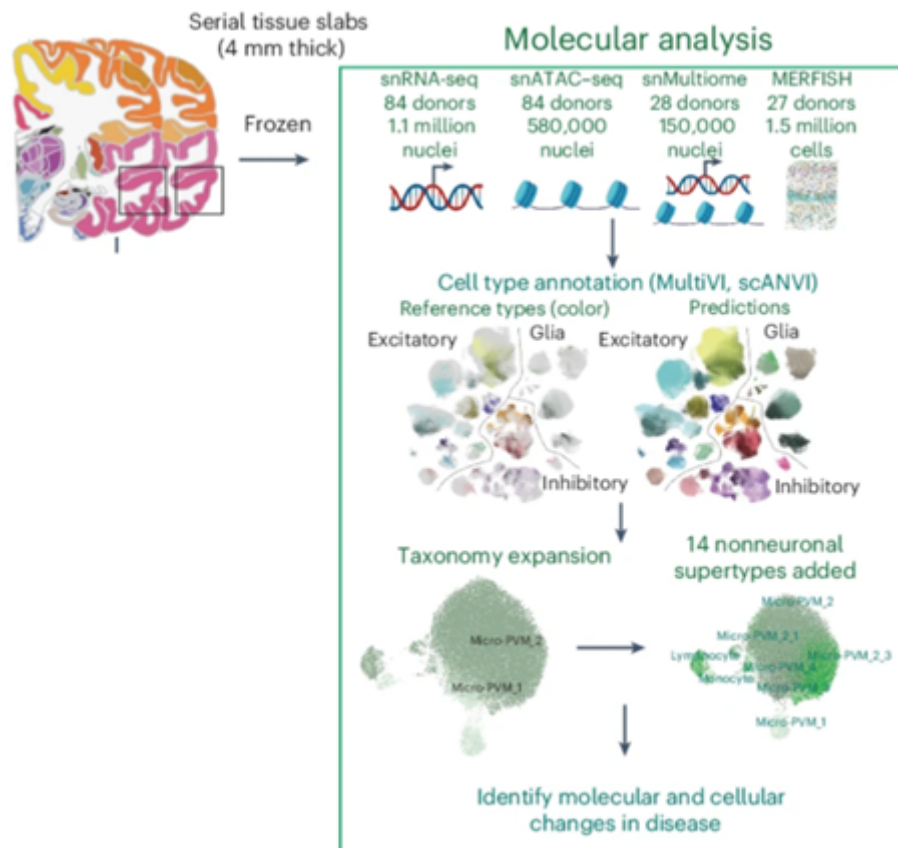


Figure 1

It might look complicated, but here's the basic idea: collect brain tissue samples, sequence RNA from

individual neurons to measure gene expression, use that transcriptional data to sort cells into different types (we'll cover this in more detail later), and then analyze the data for gene expression differences related to Alzheimer's disease

Now let's take a look at the heroic donors who were studied for the SEA-AD project. This part of Figure 1 (Panel B) gives information about the donors based on levels of symptoms for AD: high, intermediate, low, or had no AD diagnosis (the healthy controls).

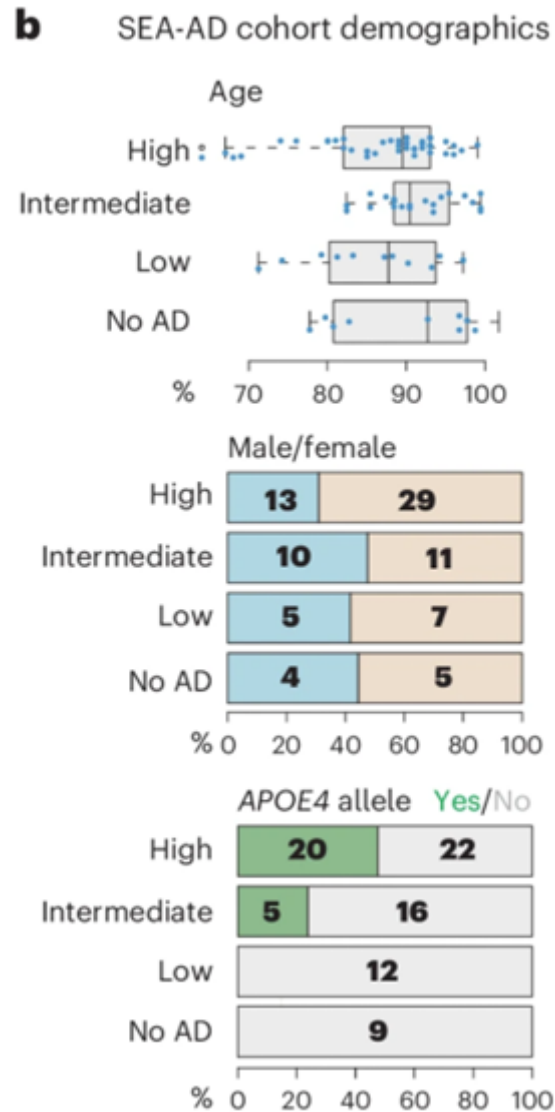


Figure 2

Check Your Understanding:

Q1: In the boxplot of ages by group above, the line in the middle of each box represents the median age. Using that median as the “typical” value, what is the typical age of donors with high levels of Alzheimer's?

1. About 70 years old
2. About 80 years old

3. About 90 years old
4. About 100 years old

Q2: What is the typical age of the donors with no AD?

1. About 70 years old
2. About 80 years old
3. About 90 years old
4. About 100 years old

Q3: Why is it important that the researchers found healthy donors of about the same age as those with high levels of Alzheimer’s disease symptoms?

1. Because if the ages were very different, it wouldn’t be possible to know if differences in gene expression were due to Alzheimer’s or to normal aging.
2. Because if the ages were very different, there could be generational differences that would make it difficult to know if the gene changes observed were really due to Alzheimer’s.
3. Because if the ages were very different, there could be other age-related factors influencing the results.
4. All of these!

Q4: Overall, were there more female or more male donors in this study?

1. More female donors
2. More male donors
3. The exact same

Q5: Female donors are only slightly more represented than male donors in the “No AD” group. Is the same true for the high AD group?

1. Yes, it’s the same proportion.
2. No, female donors are under-represented in the high AD group.
3. No, female donors are even more over-represented in the high AD group.

Understanding and Representing Gene Expression Data

The SEA-AD project measured gene expression from individual cells in the CNS of AD and healthy donors. Because RNA is transcribed from DNA this is often called transcriptomic data – data that measures all the transcripts a single cell currently expresses. Additionally, in order to classify different cell types based on gene expression, we use a representation called a UMAP (Uniform Manifold Approximation and Projection). To represent levels of gene expression of specific genes to compare across different cell types we can also use a heatmap.

Let’s hear from [Meuler & Casimo \(2024\)](#) on their explanation of how single cell RNA sequencing works and how it is represented as UMAP (Uniform Manifold Approximation and Projection) and heatmaps.

The remainder of the text and figures in this section are adapted from their work:

“Transcriptomic data is a type of data that allows scientists to investigate which genes a cell is transcribing/expressing and in what quantities. If a cell, and more specifically, that cell’s nucleus, contains a specific RNA transcript, this indicates that the cell is expressing the specific gene associated with that RNA. By (1) isolating nuclei, (2) sequencing the mRNA transcripts found within the nuclei, and (3) counting those transcripts, we can tell which genes the cell is expressing and how much these cells are expressing these genes. Repeating this process for thousands of cells from a sample of brain tissue allows researchers to find

similarities and dissimilarities between cells on the basis of their gene expression. These patterns of similarity and dissimilarity are then used to classify certain cells as specific “types.” The graphic below explains in detail how scientists gather and interpret transcriptomic data.”

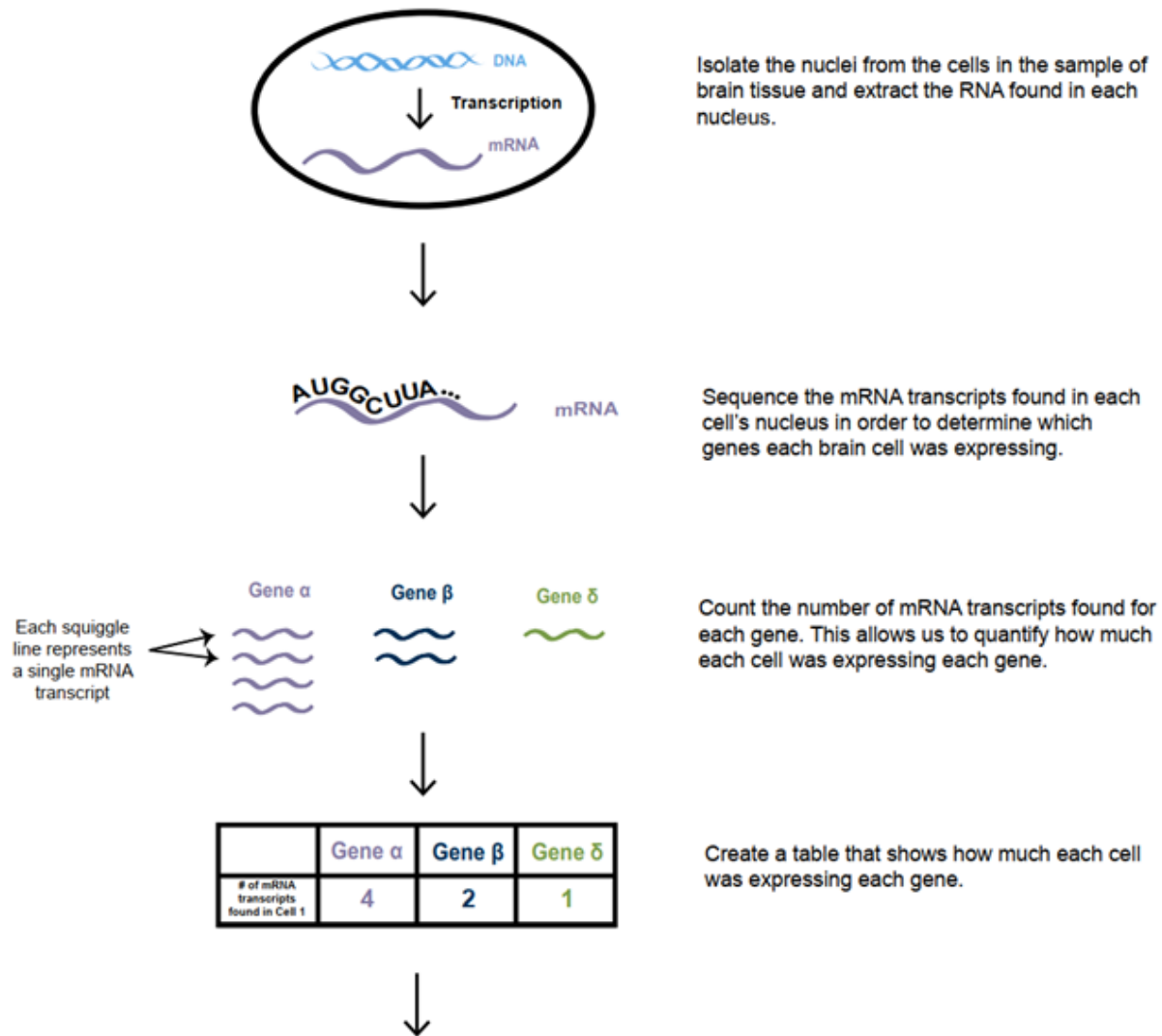


Figure 3

Check Your Understanding:

Q6: The SEA-AD project measured _____ in individual cells from the brains of healthy and AD donors.

1. DNA
2. RNA
3. Protein
1. All of the above


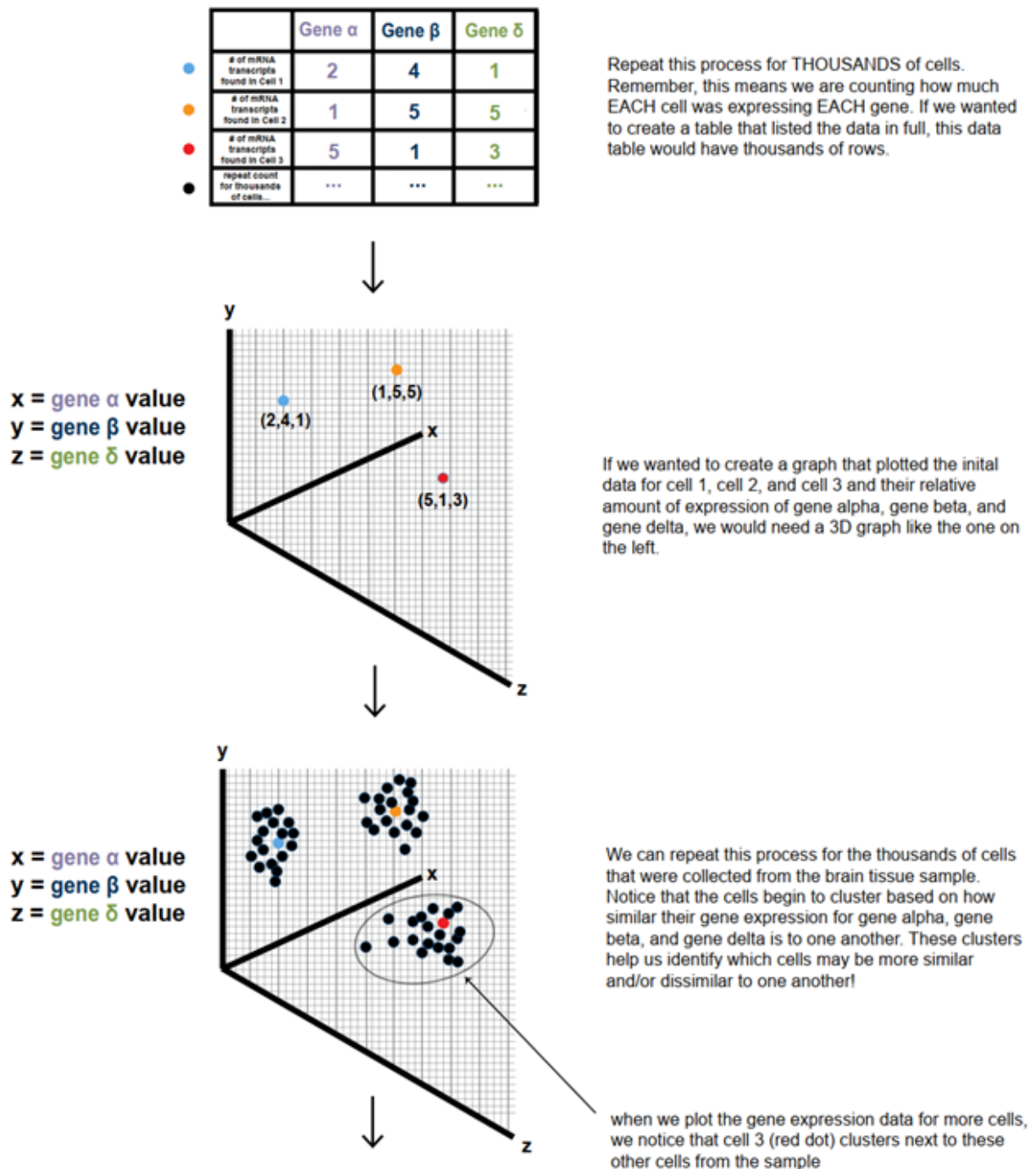


Figure 4

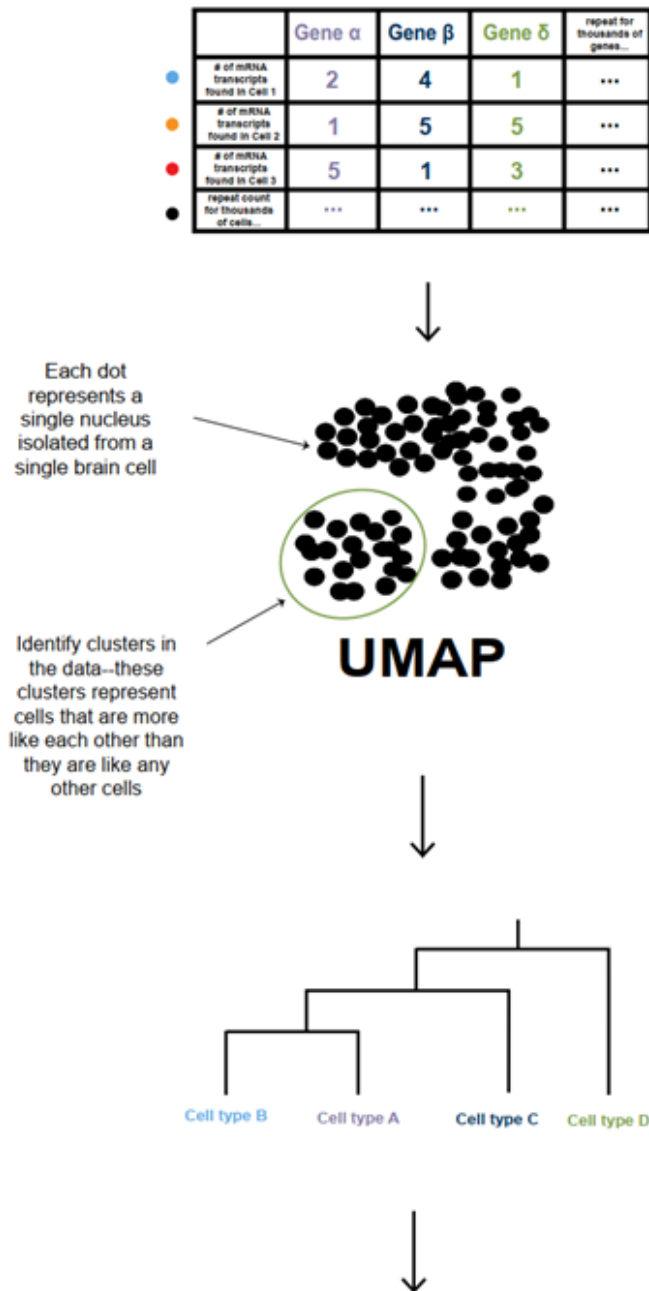

In addition to collecting data on gene expression for thousands of cells, scientists will add another layer of complexity by measuring the gene expression of these thousands of cells for THOUSANDS of genes. A table displaying this data would have thousands of rows and thousands of columns. Since the graph would now have much more than just 3 dimensions, we will need a special type of tool to graphically represent this data in a way that humans can visualize.

In order to plot this many-dimensional graph in a way humans can visualize, we use a dimensionality reduction tool, such as a UMAP, to plot it in a 2D space. Dimensionality reduction is a technique that helps represent many-dimensional data in just two or three dimensions.

Organize the clusters identified in the UMAP to construct a dendrogram that displays hierarchical relationships between the clusters based on each cell type's similarity and dissimilarity of gene expression.

Figure 5

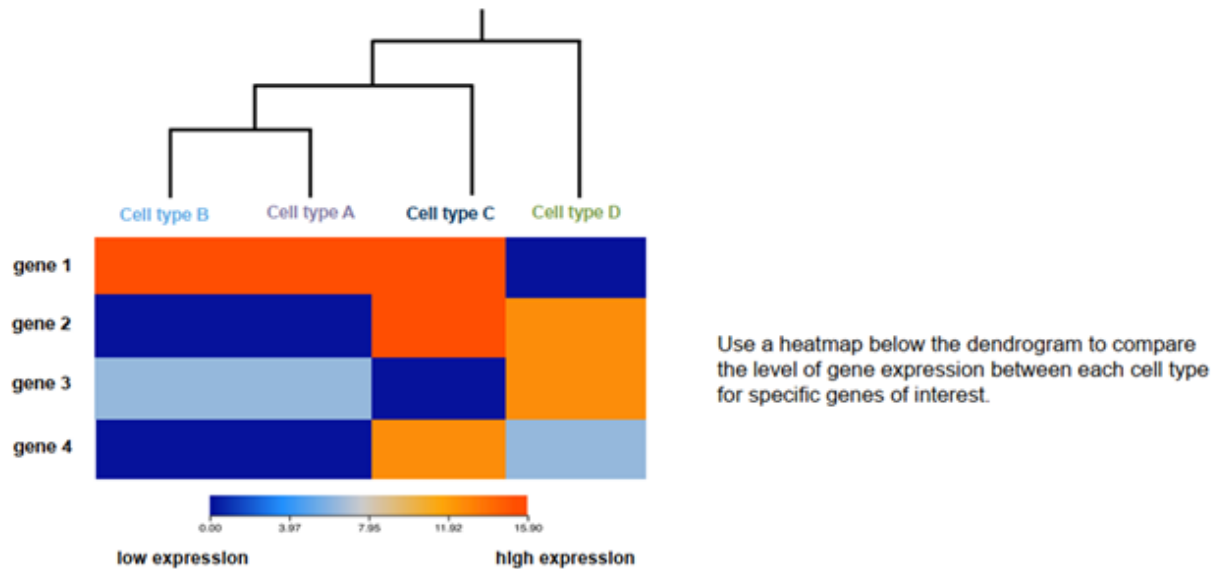


Figure 6

Q7: For each cell, sequencing enabled the SEA-AD researchers to quantify:

1. The expression of one gene for each cell
2. The expression of many genes for each cell
3. The expression of one protein for each cell
4. The expression of many proteins for each cell

Q8: For the UMAP plots shown in Figure 1 of Gabitto 2024, each dot represents a single:

1. Patient
2. Cell
3. RNA
4. Brain region

Using the Cell×Gene Explorer for the SEA-AD dataset

Let's roll up our sleeves and figure out how to explore the SEA-AD data. It's a complex data set, but we can do it!

1. **Access the SEA-AD data through the Cell x Gene explorer:** <https://cellxgene.cziscience.com/collections/1ca90a2d-2943-483d-b678-b809bf464c30>

You'll find many options to explore. The first item on the list is the complete SEA-AD dataset for the middle temporal gyrus ("Whole Taxonomy - MTG: Seattle Alzheimer's Disease Atlas (SEA-AD)"). Next, you can access the full dataset for the dorsolateral prefrontal cortex ("Whole Taxonomy - DLPFC: Seattle Alzheimer's Disease Atlas (SEA-AD)").

Following these are breakdowns of the data by specific cell types (for example, just the Lamp5 cells, just the Pax6 cells, etc.). These subsets help make the project's vast data more manageable and easier to focus on.

For our first exploration, let's start with the whole dataset for the middle temporal gyrus. Click **Explore** next to the "Whole Taxonomy - MTG: Seattle Alzheimer's Disease Atlas (SEA-AD)" option to begin.

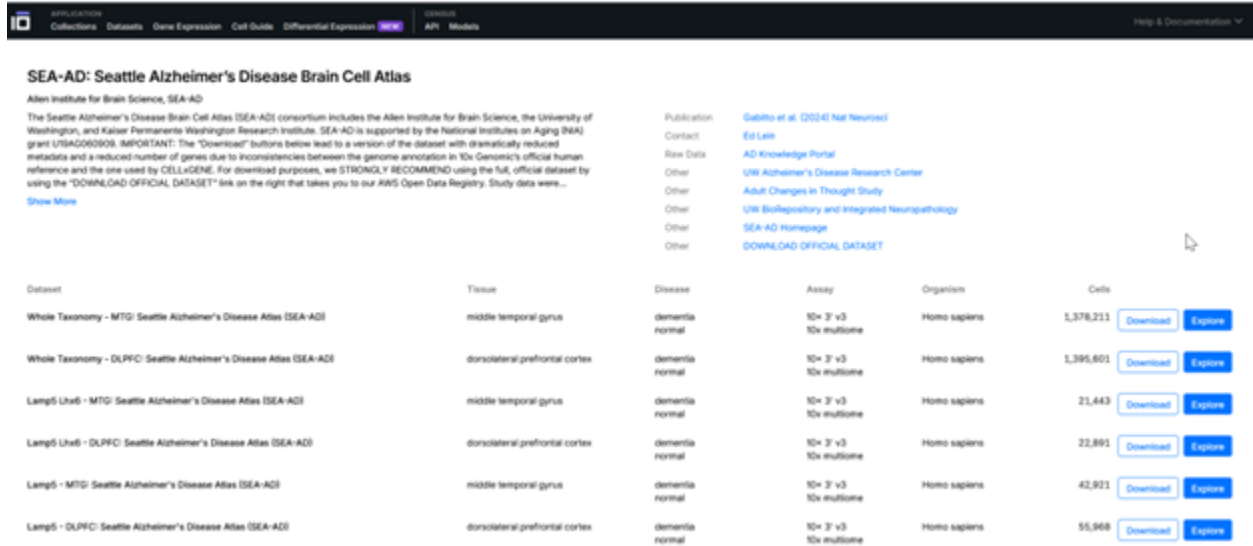


Figure 7

Dataset	Tissue	Disease	Assay	Organism	Cells
Whole Taxonomy - MTG: Seattle Alzheimer's Disease Atlas (SEA-AD)	middle temporal gyrus	dementia normal	10x 3' v3 10x multiome	Homo sapiens	1,378,211 Download Explore

Figure 8

2. **Get to know the interface.** You should now see the Cell × Gene explorer displaying the SEA-AD data. Your screen will look like this: filters and mappings on the left, a UMAP visualization of all individual cells in the center, and the gene explorer panel on the right.
3. **Use the zoom tool** - In the large middle panel, you'll see the UMAP — a visualization representing every cell sequenced in this database, totaling 1,378,211 cells! Right now, it might just look like a big black blur. Let's zoom in to see more detail! Click the Zoom tool, then use your mouse to zoom in and out of the UMAP.

As you zoom in, you'll begin to see individual small dots – each dot represents a single cell – it could be a neuron, or nonneural brain cell type such as a microglia, or an astrocyte... whatever was sequenced

Now Zoom back out to where you can see the whole UMAP.

4. **Color cells by a category with the raindrop tool:** To better understand the UMAP, we can color-code the cells based on their characteristics using the raindrop tool. Start by clicking the raindrop tool and selecting **disease**. This will color the cells so that those from donors with Alzheimer's disease appear purple, while cells from healthy donors appear green. Pretty cool, right?

Try coloring by other characteristics. For example, can you color the cells by the sex of the donor (males/females)?

5. **Understand the UMAP:** This looks cool, but why are the cells shown as these lumpy smudges? The position of each cell is determined by how similar it is to other cells based on gene expression. In other words, cells placed close together have similar gene expression profiles across their entire transcriptomes, while cells farther apart are less similar.

Placing cells on a UMAP like this reveals clusters of transcriptionally similar cells, which is how researchers in this project defined different cell types.

To explore this yourself, use the raindrop tool to color by **cell_type**. Then, expand the **cell_type** heading

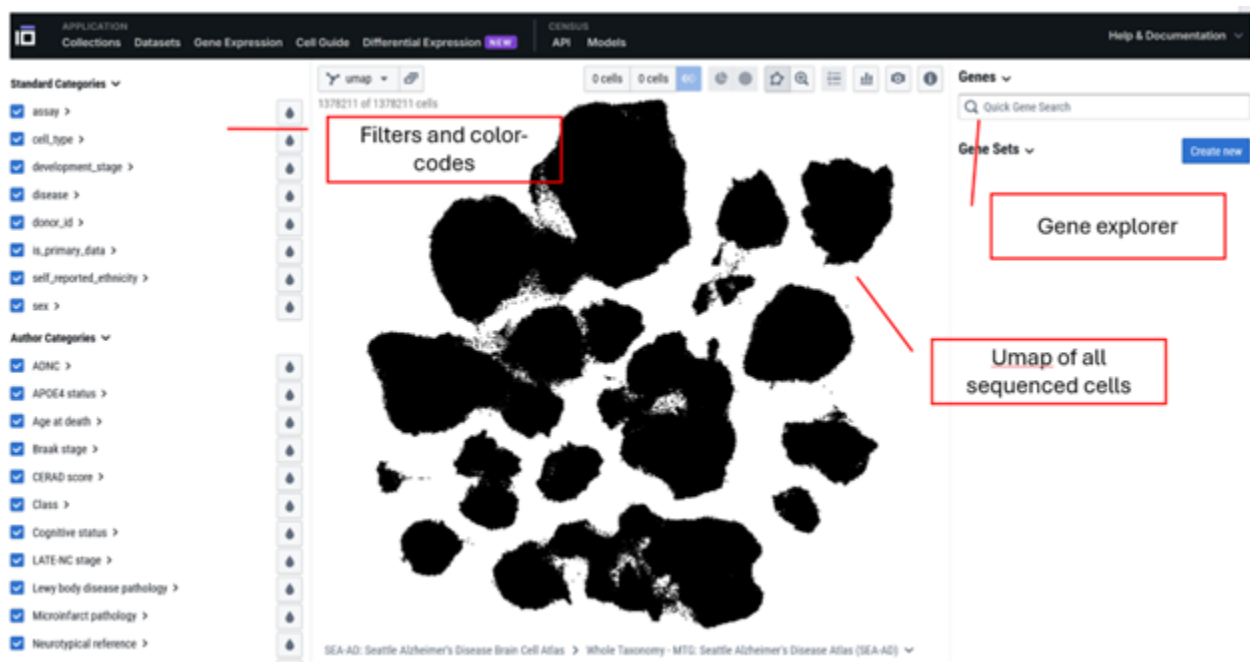


Figure 9

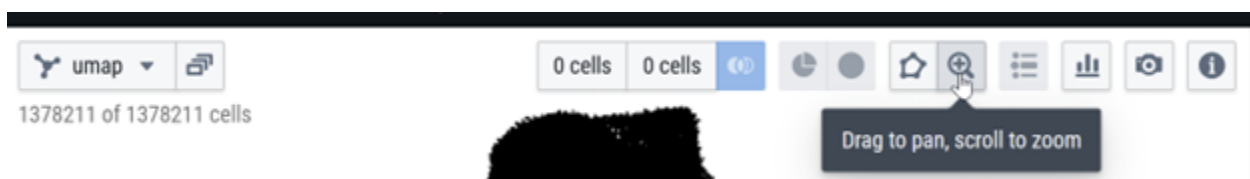


Figure 10

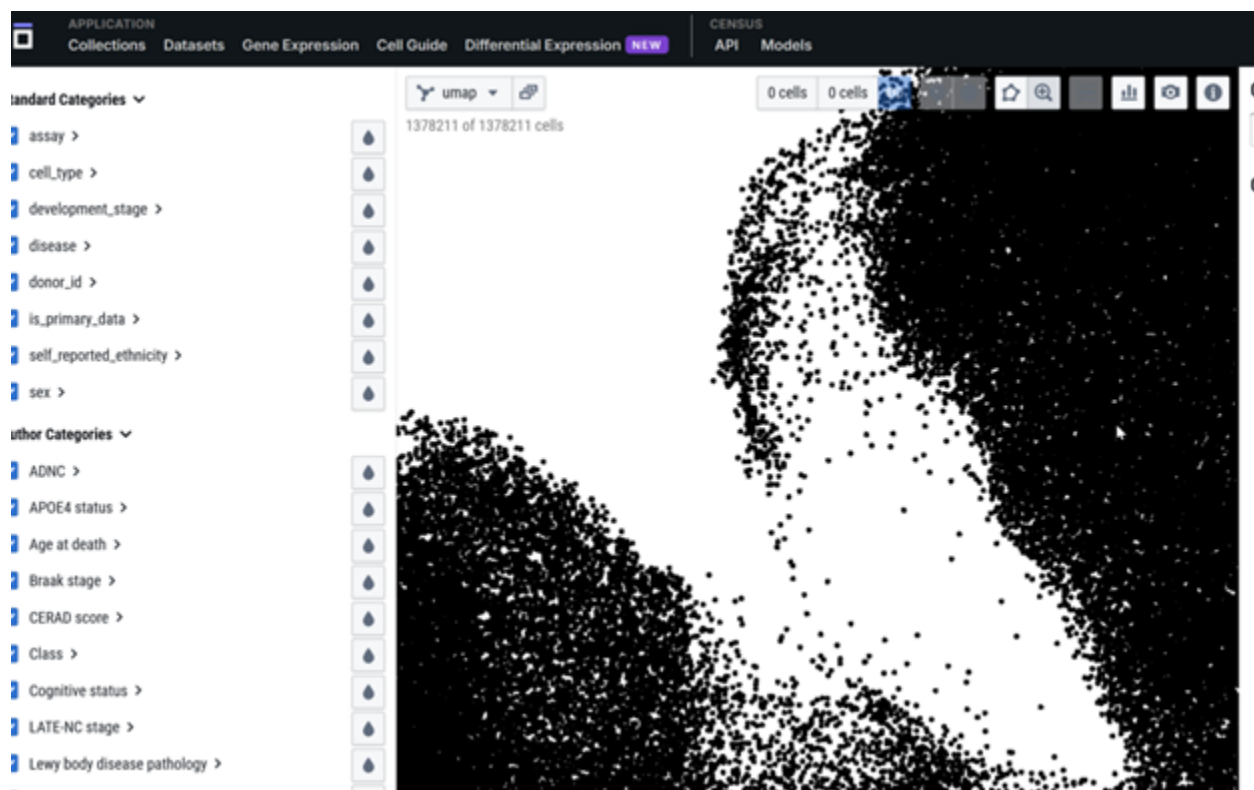


Figure 11

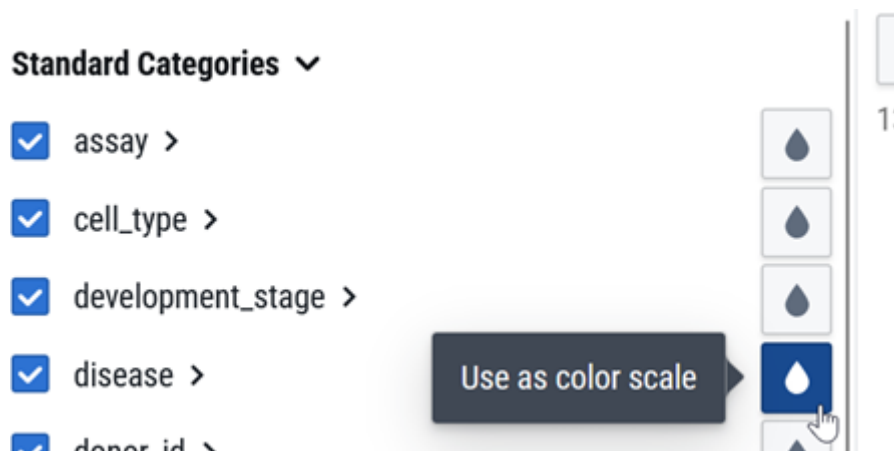


Figure 12

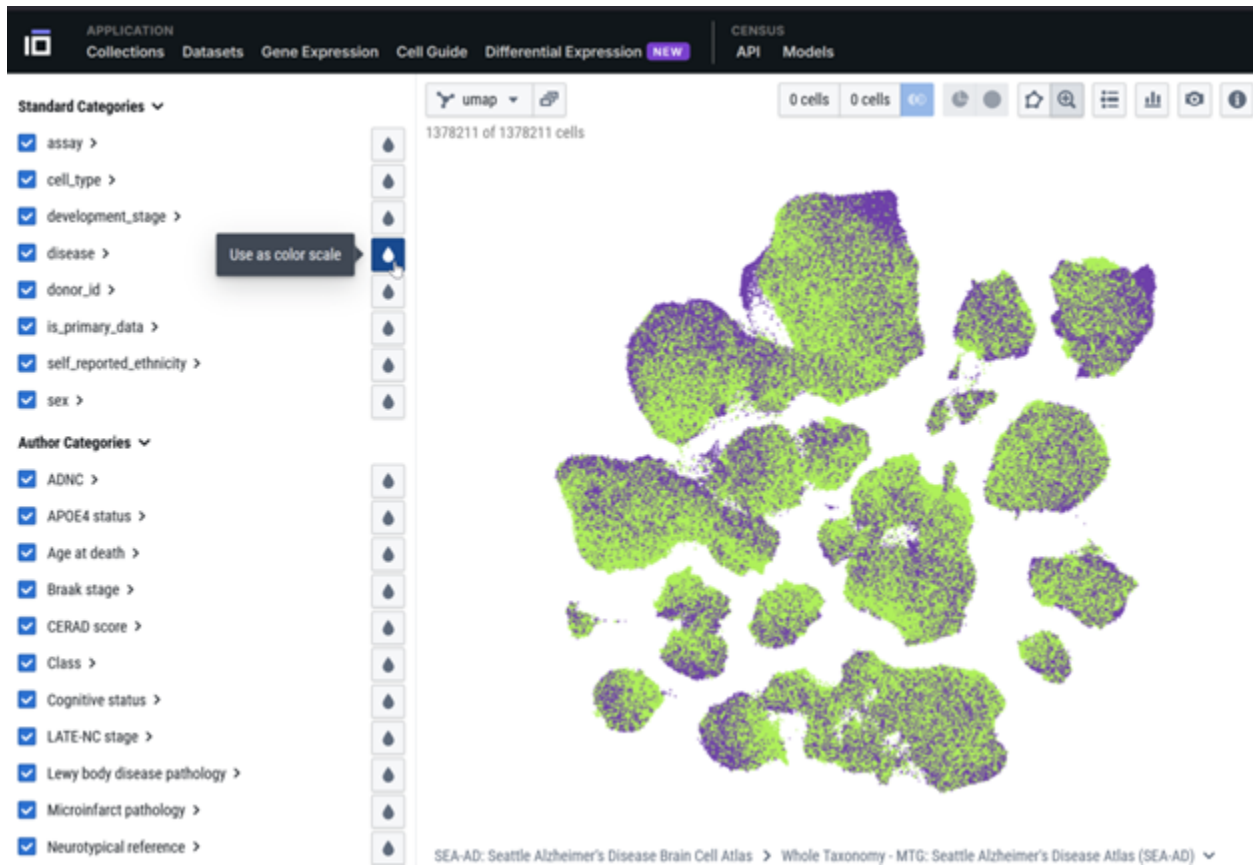


Figure 13



Figure 14

(click on the arrow >) to see a detailed view — your screen should look like this:

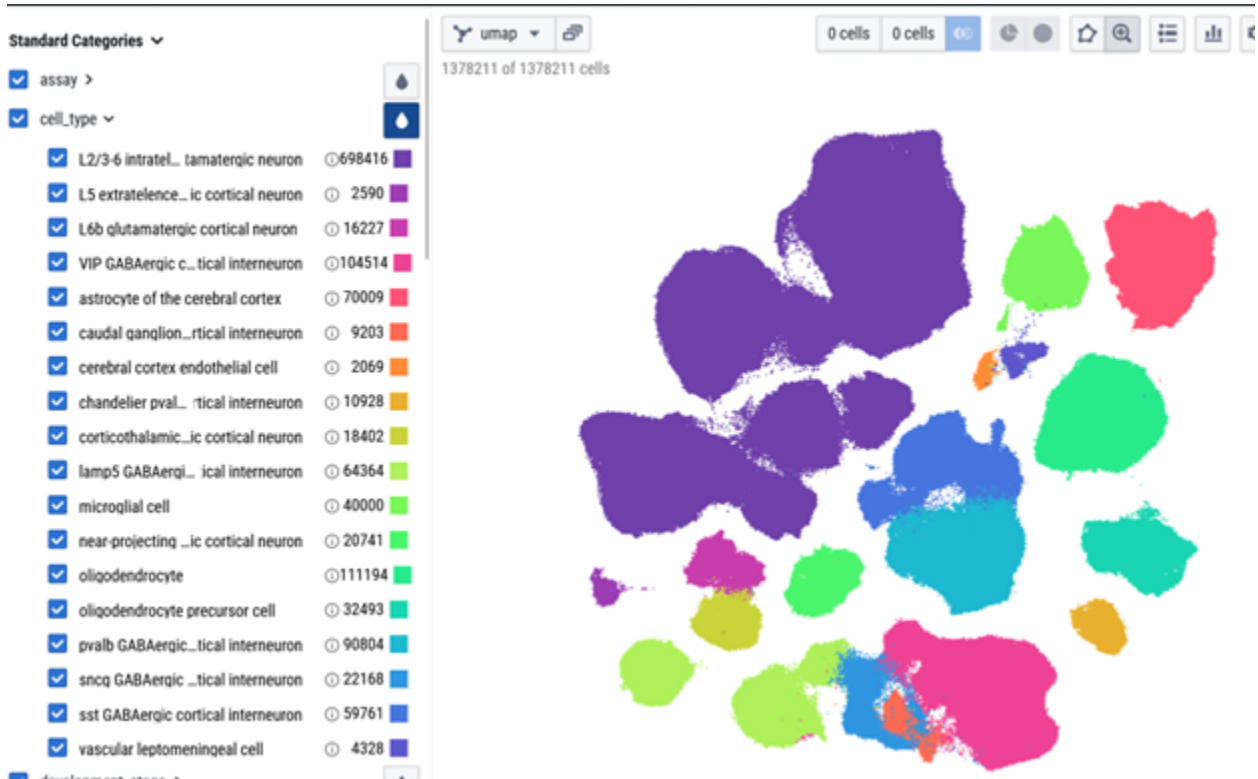


Figure 15

When you hover your mouse over different cell types in the list, they'll be highlighted on the map. To label the different populations, click the list icon located to the right of the zoom tool.

You'll notice the purple “islands” represent L2/3-6 intratelencephalic projecting glutamatergic neurons. The bright green islands at the bottom left are Lamp5 GABAergic cortical interneurons. Notice how the glutamate islands and GABA islands are quite far apart — excitatory and inhibitory neurons have distinctly different gene expression profiles!

6. **Understand the group-level summaries.** Once you've color-coded your cells, the CellxGene explorer has a neat feature — it shows how that color coding breaks down across other groups. It might sound a bit tricky, but stick with me!

Collapse the **cell_type** list and then expand the **disease** list. You'll see a stacked bar graph showing the distribution of cells in the dementia group versus the normal group.

This graph shows the breakdown of different cell types — like the purple glutamatergic neurons and the green GABAergic neurons — within the dementia group and the normal group. While the dementia group has 536,608 cells and the normal group has 841,403 cells, the stacked bars let you compare the proportions of each cell type in both groups.

Although we can't draw definitive conclusions from this summary alone, generally, what does this suggest?

- (1) That there are radical differences in the cell types found in dementia (AD) donors compared to normal donors.

OR

- (2) That both groups appear to have roughly the same breakdown of cell types.

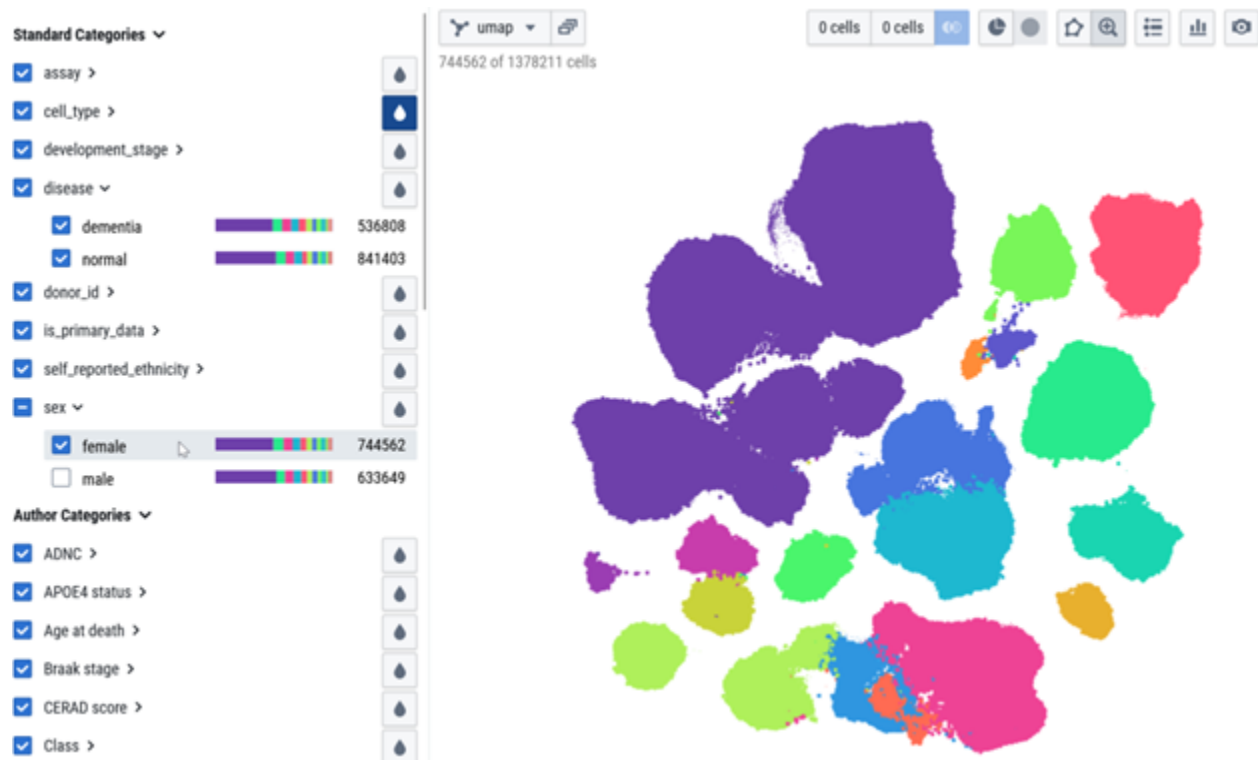


Figure 16

If you picked the second option, great job! The cell type distributions do look fairly similar, which suggests Alzheimer's may not cause dramatic changes in cell types within the middle temporal lobe. Interesting, right?

7. **Learn how to filter.** Instead of color-coding, you can also filter the data to focus on specific cells. For example, you might want to explore gene expression changes only in female donors. Filtering is simple — just expand the category you want and select the cells to keep.

Filtering lets you narrow your focus to specific groups of interest, making it easier to detect patterns or differences that might be hidden when looking at the entire dataset. This targeted approach can reveal insights specific to certain populations, like female donors, that could be important for understanding disease mechanisms. Here's an example showing only the cells from female donors:



8. **Let's start working with some genes!** If you've been filtering, be sure to reset and select all cells again. Now, we're ready to dive into gene expression.

In the Quick Gene Search box (right side of screen), type XIST and hit enter.

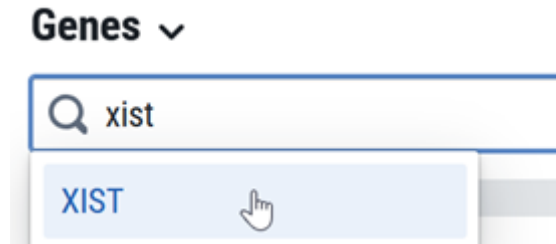


Figure 17

And add it to your gene list:

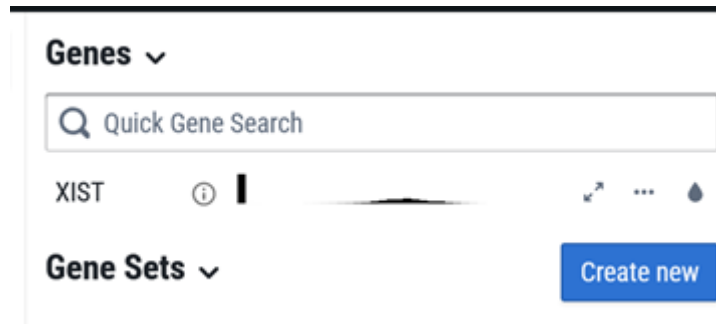


Figure 18

9. **Color-code by Gene Expression.** The XIST gene has our old friend the water-drop tool. Click it and you can color code the cells in our UMAP by their level of XIST expression.



Figure 19

Your UMAP will now look like this:

What does this show us? The green cells have very low expression of **XIST** – while they have two copies of the gene in their DNA, they aren't actively transcribing much **XIST** mRNA. In contrast, the blue cells show higher levels of **XIST** mRNA expression.

10. **Look at the overall expression histogram.** Click the <> (double arrow) button to expand the information about XIST. You'll see a histogram showing the distribution of expression levels for XIST among all the cells in the current database.

A histogram shows how common different scores are. In this case, the Y-axis represents the number of cells, and the X-axis shows gene expression levels. Each colored bar tells you how many cells have that specific level of expression for the gene.

Think of it like a popularity chart at a party – the taller the bar, the more people (or cells) have that particular “popularity score” (gene expression level).

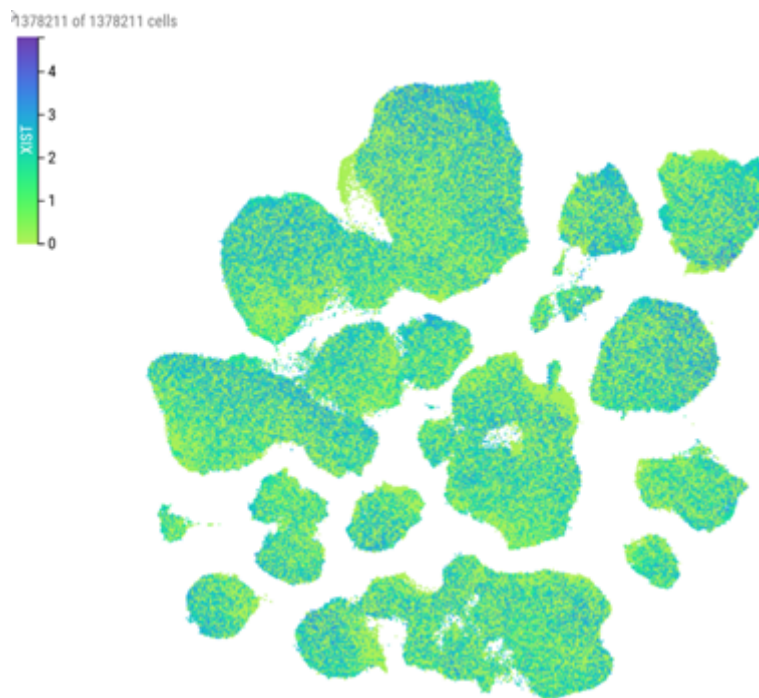


Figure 20



Figure 21

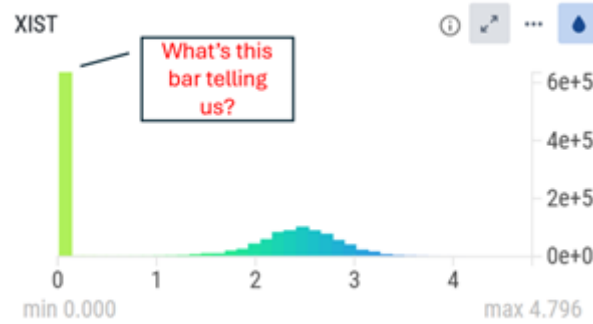


Figure 22

This histogram looks unusual! Let's try to make sense of it.

First, notice the huge bar right at zero. What does that tell you?

- That there are zero cells that fail to express **XIST**.
OR
- That there are a large number of cells that don't express **XIST** at all.

If you chose the second option, great job! That big bar shows there's a large group of cells with almost no **XIST** expression.

Now, take a look at the other cluster of bars forming a nice "bell curve":

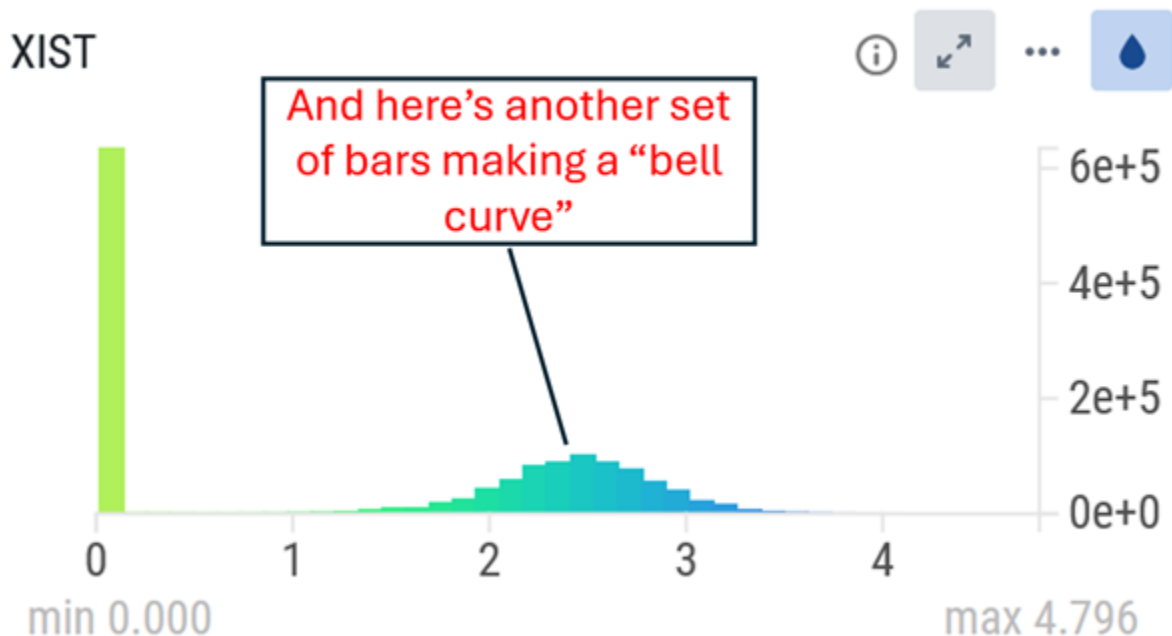


Figure 23

What does this other set of bars tell us?

- That there is a group of cells expressing **XIST**, but at varying levels—some a little, some more.

OR

- That **XIST** must encode a protein shaped like a bell.

This histogram reveals two distinct patterns of **XIST** expression: a large group of cells that don't express it at all, and another group that does, with varying amounts.

Q9. Why do you think there's this difference? What distinguishes these two groups of cells?

11. **Find out more about a gene.** Click the (i) icon to find out more about **XIST**. You'll get a brief summary as well as a link to NCBI's information about this gene.

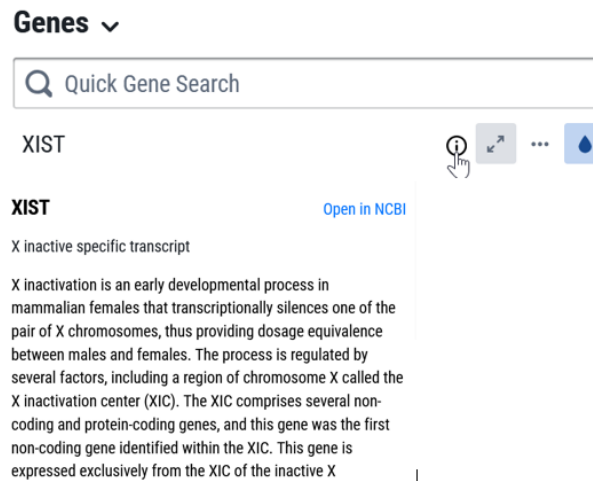


Figure 24

Q10. From this description you might hypothesize that:

- A. **XIST** is a gene only expressed in females
- B. **XIST** is a gene only expressed in males
- C. **XIST** is a gene of unknown function, so maybe it is the cause of AD

What this tells us is that **XIST** is a gene responsible for silencing one of the two X chromosomes in genetic females. Ah-ha! This means **XIST** is expressed in female cells but not in male cells – which helps explain our histogram, where some cells show no **XIST** mRNA and others show quite a bit.

Let's dig in and see if that's really the case!

12. Explore Differential Expression by Comparing Expression Maps Across Groups

To learn more about **XIST** expression, we'll compare its levels in cells from female donors versus male donors. The CellxGene explorer isn't the simplest tool for this, but don't worry — we will work through it together.

First, roll down the selector for sex on the filter/mapping panel on the left:

Right away, the CellxGene explorer is showing you the **XIST** expression broken down by sex.

Q11. What can you tell from this?

- A. That cells from female donors are the ones expressing **XIST**, corresponding to the “bell curve” in the histogram of all cells.
- B. That cells from male donors do not express **XIST**, matching the large peak at zero in the histogram.
- C. That **XIST** expression is highly sex-dependent: cells from male donors show little to no expression, while female donor cells express **XIST** at varying levels.



Figure 25

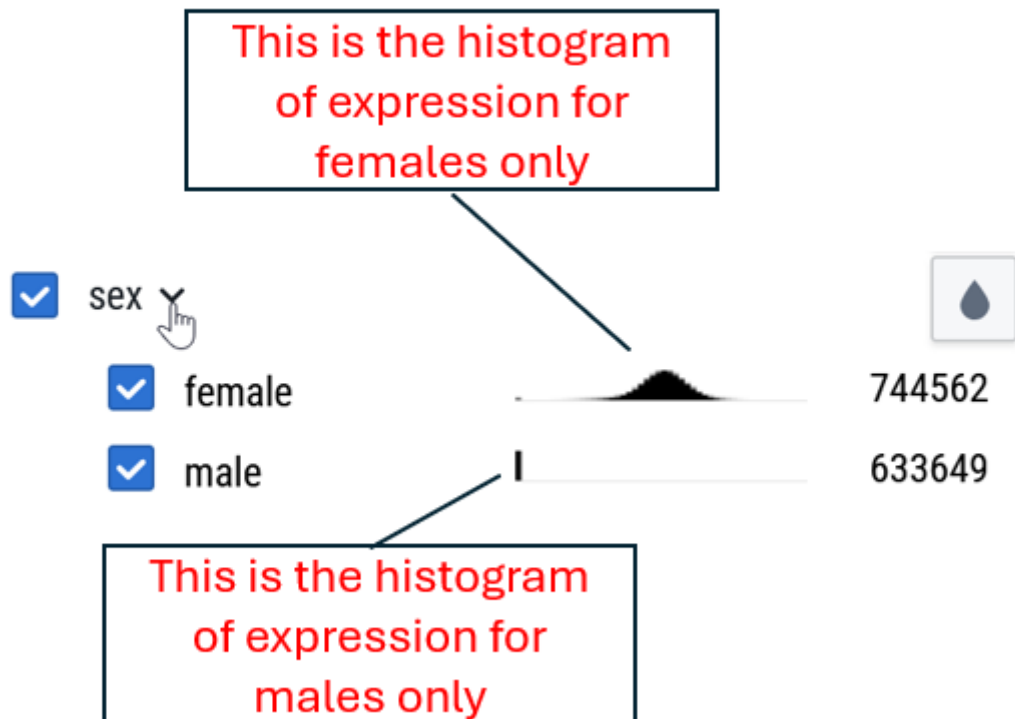


Figure 26

D. All of the above.

You can see that female cells display that nice bell-shaped curve of expression—most produce **XIST** mRNA, with some expressing a bit more or less. In contrast, male cells only show the tall bar at zero, indicating almost no **XIST** expression.

This is a qualitative difference in gene expression – if you picked a single cell and checked its **XIST** level, you’d have a pretty good idea whether it came from a biological female or male donor.

Wouldn’t it be great if we could find such clear-cut gene expression differences related to Alzheimer’s disease? Let’s find out as we continue the lab!

13. Make an Image Demonstrating Differential Expression

To document our findings, we’re going to create a figure that visualizes the difference in **XIST** expression between cells from male and female donors.

Here’s the plan:

1. Filter the data to show **female-only** cells (the UMAP should take on a blueish hue indicating moderate **XIST** expression).
2. After filtering, hover your mouse over **female** in the categories list to ensure it’s selected and not highlighting other groups.
3. This view shows **XIST** mRNA expression in cells from female donors.
4. You can capture this image by taking a screenshot—or better yet, use the **SNAPSHOT** tool within the explorer to save a high-quality image.

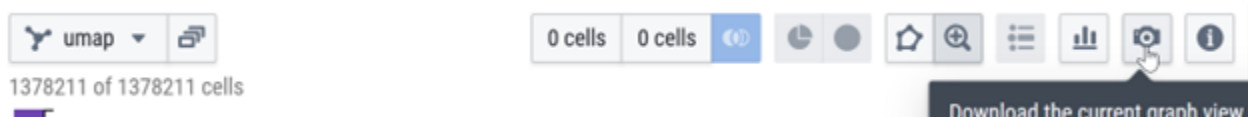


Figure 27

1. This will likely save to your Downloads folder or wherever you have set this on your computer/device
2. Note that the file names that are saved with generic names: CELLxGENE_umap_emb.png
3. Be sure to rename the saved files so you can identify which group is which!
5. Now, switch the filter to show **male-only** cells (the UMAP should turn bright green, indicating low **XIST** expression!). Again, hover your mouse over **male** in the category list to ensure no other groups are highlighted. Then, take a snapshot of this view.
6. Next, bring your two images together – one for females and one for males – and **label them clearly**. You can use PowerPoint, Google Slides, or any other tool you prefer.
7. Add a scale bar to your figure: Since the CellxGene explorer’s snapshot doesn’t include a scale bar (which isn’t ideal for scientific figures), be sure to also take a screenshot of the scale separately and include it in your combined image.

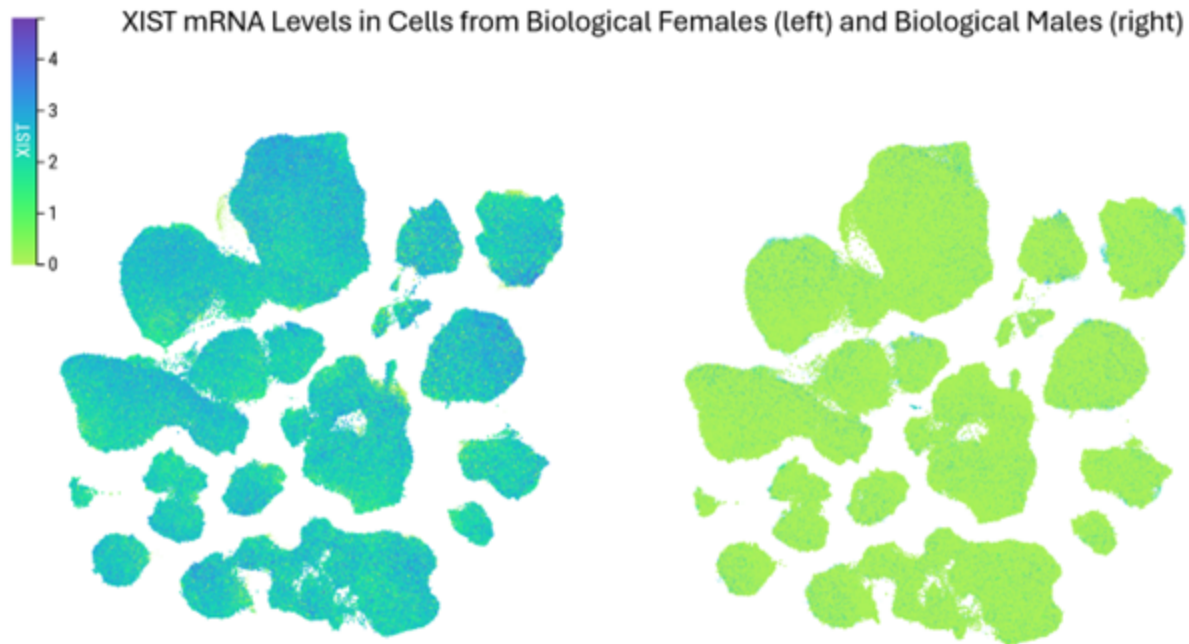


Figure 28

Your overall image should end up looking like this:

Q.12 Submit your final figure

Q.13 What can you conclude from this (your figure)?

1. That there is a qualitative difference in XIST mRNA in cells from the brains of males and females – cells from females make XIST mRNA (to varying extents), male cells largely do not (without much variability).
2. That the differential expression of XIST seems consistent across all the cell types measured in the Middle Temporal Gyrus – it is the same pattern in the glutamate clusters, the GABA cluster, etc. etc.
3. All of the above

Nice work! You've gotten hands-on experience using the Cell×Gene explorer to dive into the SEA-AD database. Below are a few challenges to test your skills in navigating and interpreting this rich dataset.

Try On Your Own

Now that you have a good handle on the tool, explore the expression of **ZNF536**. Examine **ZNF536** expression across different cell types by expanding the **cell_type** category and inspecting the histogram breakdowns. Try checking and unchecking various cell types.

Think about the following questions as you are checking/unchecking various cell types:

- Are there some cell types that express **ZNF536** strongly?
- Are there others with little or no expression?
- Can you identify any patterns that explain which cell types express this gene?

Q.15: Create a figure comparing expression between two cell types: one with high ZNF536 expression and one with low or no expression. Combine your images into a single, well-labeled slide and submit the figure.

Q.16: Describe what you can conclude from this figure in a few sentences.

You've completed the pre-lab – great job! Below are some extra resources that might be helpful. See you in lab!

More Resources

- Feeling a bit shaky about Cell×Gene explorer? There is an even more extensive tutorial available here: https://cellxgene.cziscience.com/docs/04___Analyze%20Public%20Data/4_1___Hosted%20Tutorials
- Want to read the paper by Gabitto et al. (2024) describing the SEA-AD project? It is an open-access article (hooray!) available here: <https://www.nature.com/articles/s41593-024-01774-5>

Here is a video from the Allen Institute on the SEA-AD project: <https://www.youtube.com/watch?v=lcS3o1N22Bc&t=2s>

Gabitto, Mariano I., Kyle J. Travaglini, Victoria M. Rachleff, Eitan S. Kaplan, Brian Long, Jeanelle Ariza, Yi Ding, et al. 2024. “Integrated Multimodal Cell Atlas of Alzheimer’s Disease.” *Nature Neuroscience* 27 (12): 2366–83. <https://doi.org/10.1038/s41593-024-01774-5>.